

WE ARE ALL UNIQUE.

Kodak Unique HD Progressive Lens

WE ARE ALL UNIQUE.



New Levels of Customization

When first introduced, **Kodak** Unique[™] Progressive Lens quickly became an established favorite lens design. With its ease of dispensing and high rate of adaptation, **Kodak** Unique Lens truly lives up to the expectation of being *unique*.

With this proven success, we are now able to introduce the next level of **Kodak** Unique Lens customization with **Kodak** Unique HD Lens. The HD stands for 'High-Definition' and fulfills this description by including additional measurements to more highly customize the lens to the individual patient's viewing needs.

We are all unique so why shouldn't our viewing experience be as well?




See Life in High Definition

Kodak Unique HD Lens includes the features of Kodak Unique Lens:

- ▶ Full backside progressive design
- ▶ Vision First Design[™] Technology
- ▶ i-Sync™ Technology
- Six corridor lengths
- Variable decentration
- Wide distance vision
- Prescribed prism

+ more...

Creating a Unique Wearing Experience

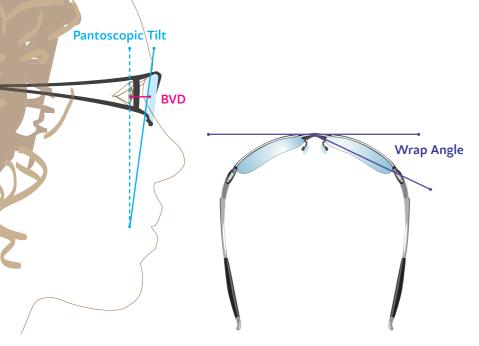


Kodak Unique HD Lens customization includes the following additional features:

Prescription Compensation

The patient's **Back Vertex Distance (BVD)**, **Pantoscopic Tilt** and **Wrap Angle** are included to provide a customized lens optimized for the individual patient's prescription and frame selection.

Variable Inset

Relying on the Near Reading Distance, Lens Power, Back Vertex Distance (BVD) and Monocular PD (MPD), Kodak Unique HD Lens delivers a near viewing experience tailored to the individual patient's needs.

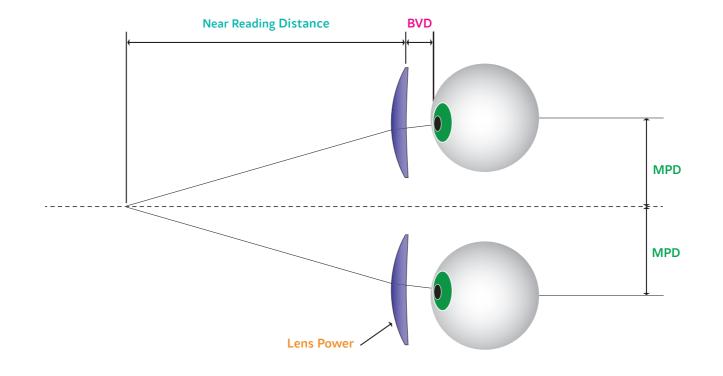

The Patient's Viewing Experience

- Wide fields of view offer picture-perfect vision at all distances
- Advanced technologies allow greater visual comfort
- Reduced distortion in the design for smooth, natural viewing and comfortable wear
- Greater color contrast for a more dynamic viewing experience



> Thinner lenses resulting in more cosmetically-appealing eyeglasses

The ECP's Dispensing Experience



- Outstanding Patient Adaptation
- A Trusted Brand Name
- Easy to Dispense
- Advanced Technologies at a Sensible Price
- A Variety of Material Options



Kodak LENS INTELLIGENT DISPENSING SOFTWARE

Vision First Design Technology[™] is a proprietary technology used to produce **Kodak** Unique HD Progressive Lenses.

Vision First Design Technology provides:

- A broad field of view in a clear distance area
- Smooth gradation of power across the surface of the lens to ease patient adaptation
- Gentle binocular balance for quick, clear object recognition
- Visual quality often compared to a single vision lens

i-Sync

i-Sync[™] Technology elevates the level of optical performance by maximizing clarity in the peripheral areas of the lens and reducing distortion to the lens edge.

i-Sync Technology provides:

- More consistent optical performance over the range of prescription power
- Wider near viewing areas for hyperopes
- Improved distance area for myopes
- Improved image quality in principal viewing areas
- Flatter base curve capabilities

Kodak Lens Intelligent Dispensing Software

Compatible with iPad Air and iPad with Retina display (4th generation), **Kodak** Lens Intelligent Dispensing Software (IDS) is an app composed of three separate sections: FRAME, LENS and MEASURE.

FRAME: Compares and shares photos of patient frame options.

LENS: Assists the eyecare practitioner in educating the patient about the visual benefits of **Kodak** Progressive Lenses including **Kodak** Lens designs as well as polarized, photochromic and AR coating options.

MEASURE: Optimizes the use of an iPad by quickly turning it into a measuring device – providing monocular PD measurements for far and near vision, monocular fitting height measurements, pantoscopic tilt, back vertex distance, face form angle, A, B and DBL measurements.

Check out the **Kodak** Lens IDS informational video on the Signet Armorlite YouTube channel, www.youtube.com/SignetArmorliteUS.



Feature	KODAK Unique HD Lens
Corridor Lengths	6
Material Availability	50+
Full Backside Progressive Des	ign 🗸
Vision First Design Technolog	у 🖌
i-Sync Technology	✓
Variable Decentration	✓
Wide Distance Vision	 ✓
Prescribed Prism	 ✓
Prescription Compensation	✓
Variable Inset	✓

Kodak Unique HD Progressive Lens begins with the respected **Kodak** Unique Lens design and adds a higher level of prescription customization never seen before in a **Kodak** Progressive Lens.

KODAK Unique HD Lenses Because We Are All Unique.

For over sixty-six years, Signet Armorlite has been a leader in the optical industry with quality products and innovative technologies.

Signet Armorlite, Inc. 800.950.5367 www.signetarmorlite.com



