


Shamir Duo™

A New Visual Harmony

TECHNICAL INFORMATION

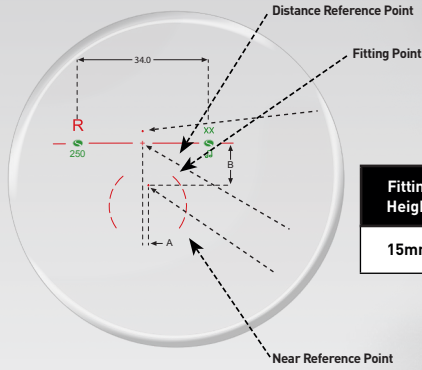

Fitting Height	"A" (see diagram)	"B" (see diagram)	"M" Invisible Marking (see diagram)
15mm	1.5	11	♪

Minimum Fitting Height: 15mm
 Ink Marking: Red • Invisible Marking: Green

Material Availability:

Material	Power Range	Add Power	Cylinder
1.5 Hard Resin (C,TV,T7,P,TX,D)	-9.00 to +6.00	0.75 to 4.00	to -6.00
DLC™ Trivex™ (C,TV,T7,P,TX,BZ)	-13.25 to +6.00	0.75 to 4.00	to -6.50
1.56 (B)	-14.00 to +6.00	0.75 to 4.00	to -6.75
NXT™ (PH)	-1.75 to +6.00	0.75 to 4.00	to -6.50
1.6 SuperLite™ (C,T7,P,TX)	-15.00 to +6.00	0.75 to 4.00	to -7.00
1.67 SuperLite™ (C,T7,P,TX,BZ)	-16.75 to +8.00	0.75 to 4.00	to -8.00
Polycarbonate (C,TV,T7,P,TX,D,BZ)	-14.75 to +6.00	0.75 to 4.00	to -7.00
1.74 SuperLite™ (C,T7)	-18.75 to +12.00	0.75 to 4.00	to -9.00

POWER RANGE IS DETERMINED BY SHAMIR CERTIFIED BLANK RANGE PRISM TO 10 DIOPTERS IN EACH QUADRANT

LEGEND:

- C-Clear
- T7-Transitions® Signature7™
- PH-Photochromic
- P-Polarized
- TV-Transitions® Vantage™
- TX-XTRActive™
- B - BluTech
- D - Drivewear
- BZ-Blue Zero™

WARNING: Polycarbonate contains a chemical known to the state of California to cause birth defects or other reproductive harm.



How to Order Shamir Duo™

Provide the following information:

- Patient prescription
(Distance correction plus Addition Power)
- Mono PD
- Fitting height *(measured from the center of the pupil to the lower rim of the frame)*
- Frame data: A, B, ED and DBL
- Frame type

Fitting Instructions:

- The Shamir Duo™ fitting point should coincide with the center of the pupil, using the invisible markings for verification (just like progressive lenses)

Technical Inquiries: 888.370.0736



[www.facebook.com/
ShamirInsight](http://www.facebook.com/ShamirInsight)



[www.twitter.com/
ShamirInsight](http://www.twitter.com/ShamirInsight)

DUO-TEC-2PA-020118



SHAMIR

ReCreating Perfect Vision®

www.shamir.com 877.514.8330