


# Shamir Spectrum+™

Everyday Lenses Enhanced for Comfortable Reading

## TECHNICAL INFORMATION

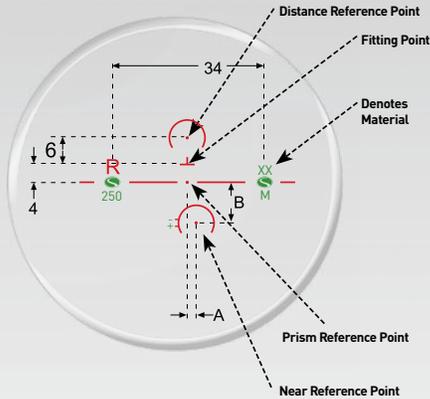

Fitting Height	"A" (see diagram)	"B" (see diagram)	"M" Invisible Marking (see diagram)
14mm	2	7.5	SP+14
16mm	2	9.5	SP+16
18mm	2.5	11	SP+18

Minimum Fitting Height:  
14mm, 16mm, 18mm  
Ink Marking: Red • Invisible Marking: Green



## Material Availability:

Material	Power Range	Add Power	Cylinder
1.5 Hard Resin (C,T,V,T7,P,TX,D)	-12.25 to +6.00	0.75 to 3.50	to -6.00
DLC™ Trivex™ (C,T7,P,TX,BZ)	-13.25 to +6.00	0.75 to 3.50	to -6.50
NXT™ (PH)	-1.75 to +6.00	0.75 to 3.50	to -6.50
1.56 (B)	-14.00 to +6.00	0.75 to 3.50	to -6.75
1.6 SuperLite™ (C,T7,P,TX)	-15.00 to +6.00	0.75 to 3.50	to -7.00
1.67 SuperLite™ (C,T7,P,TX,BZ)	-16.75 to +8.00	0.75 to 3.50	to -8.00
Polycarbonate (C,T,V,T7,P,TX,D,BZ)	-14.75 to +6.00	0.75 to 3.50	to -7.00
1.74 SuperLite™ (C,T7)	-18.75 to +11.50	0.75 to 3.50	to -9.00

POWER RANGE IS DETERMINED BY SHAMIR CERTIFIED BLANK RANGE  
PRISM TO 10 DIOPTERS IN EACH QUADRANT

### LEGEND:

- C-Clear
- T7-Transitions® Signature7™
- PH-Photochromic
- P-Polarized
- BZ-Blue Zero™
- TV-Transitions® Vantage™
- TX-XTRActive™
- B-BluTech (Indoor & Outdoor)
- D-Drivewear

**WARNING:** Polycarbonate contains a chemical known to the state of California to cause birth defects or other reproductive harm.



## Advanced Technology



### NATURAL POSTURE™ New Ergonomic Design Concept



The near viewing zone in Shamir **Spectrum+**™ is dynamically located within the lens corridor, taking into account the plus or minus power of the prescription. This reduces the need to tilt the head forwards or backwards or alternatively to lower or raise a book when reading to find the right viewing distance and angle. Any discomfort within their posture is thereby minimized and patients are ensured a comfortably natural posture during near viewing.

## How to Order Shamir **Spectrum+**™

### Provide the following information:

- Supply the Doctors refracted Rx
- Frame Data: (A, B, DBL) or frame tracer file
- Fitting height & Mono PD (far)

**Min Fitting Heights:** 14mm, 16mm, 18mm

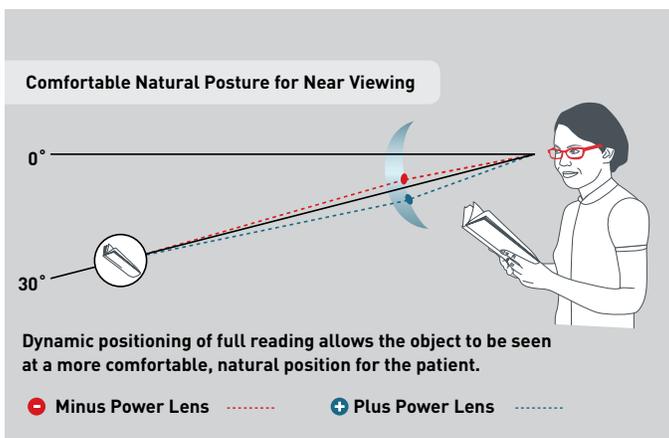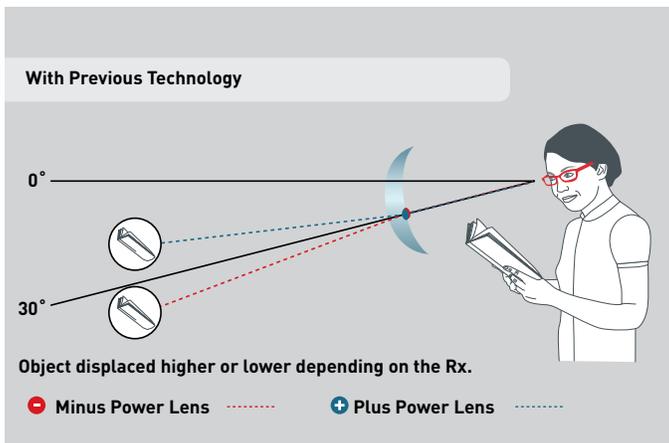

Technical Inquiries: 888.370.0736



[www.facebook.com/ShamirInsight](http://www.facebook.com/ShamirInsight)



[www.twitter.com/ShamirInsight](http://www.twitter.com/ShamirInsight)

SP+TEC-2PA-1020118



SHAMIR

ReCreating Perfect Vision®

[www.shamir.com](http://www.shamir.com) 877.514.8330